



08186 Lliçà d'Amunt  
Barcelona  
Spain

Tel.:+ 34 93 860 90 00  
Fax:+ 34 93 860 90 17  
e-mail: biokit@biokit.com  
www.biokit.com

**Date:** 27/04/05


**Ref:** 190/MKT/Q/67

**SUBJECT** **BQA – Hardware-waste separation lines**

We would like to inform you that a new BQA hardware modification has been performed to separate ISE waste solution from the rest of the waste lines.

- Quick connector for high concentrated waste has been replaced for a plastic connector.
- ISE waste line is a separated line which includes an extra labelled silicone tube and a plastic connector.

See attached picture for further details.

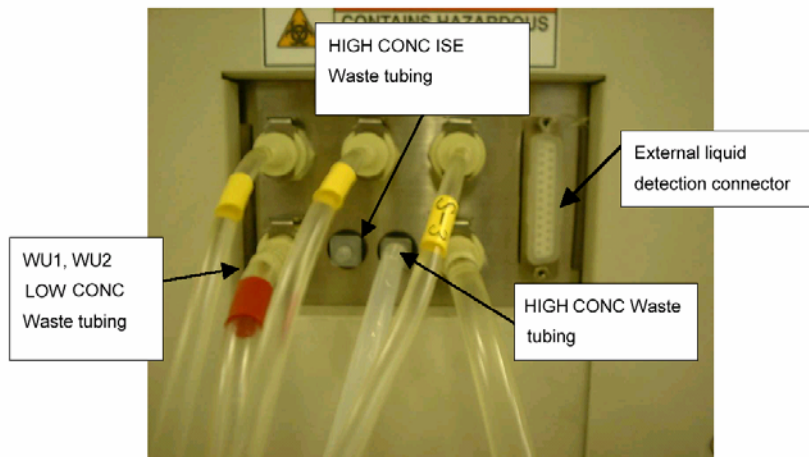

Please, don't hesitate to contact us if you need more information.

